

# Medical Pendants - Multi-movement



Precision UK supplies a comprehensive range of high quality ceiling suspension systems which have been developed under the following criteria:

- Increased Safety
- Increased Hygiene and Sterility
- Ease of use


The many advantages in the field of safety and hygiene are;

- No more trailing leads or hoses to restrict movement within the department, all equipment is attached via short connecting hoses or cables to integral hook up points.
- Fully encapsulated system, incorporating modular service boxes, to ensure maintenance and general use is unproblematic.
- Single or two piece ondal arm assemblies, motorised or static, equipment can now be located at point of use with simplicity.
- Fully equipped as desired, with high and low voltage supplies, medical gas outlet points, communication points, miscellaneous.
- Factory tested supplied to customer in modular format **certified and compliant CE 0086 Medical Device Class IIb**

The Precision UK Pendant is a modular unit, which can be supplied with or without height adjustment. The pendant consists of a ceiling mount, upon which a wide selection of arms can be mounted.

## Designed to meet the requirements of:

The UK DHSS Health Technical memorandum No. 02-01. The standard medical gas terminals installed in the pendant conform to HTM No. 02-01 and to EN ISO 9170-1. The pendant can be constructed to other standards, and other makes of terminal can be fitted.

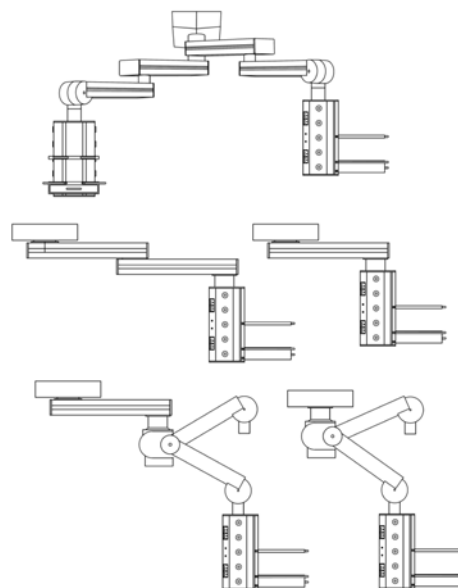

Product Code	Multi Movement Pendant Type
PEND-MM-SG	Single Arm
PEND-MM-VT-SG	Single Arm with Vertical Movement
PEND-MM-DB	Double Arm
PEND-MM-VT-DB	Double Arm with Vertical Movement
PEND-MM-TW-DB	Twin Arm, Double Swivel

For more information visit: [www.precisionuk.co.uk](http://www.precisionuk.co.uk)